


# Wheelchair scale with handrail and RS-232 interface

## seca 677 r

- + Large wheelchair accessible platform
- + Convenient and foldable to save space
- + Mobile due to transport wheels
- + Handrail as a standing aid
- + Non-slip rubber mat
- + With RS-232 interface for transmitting weighing results

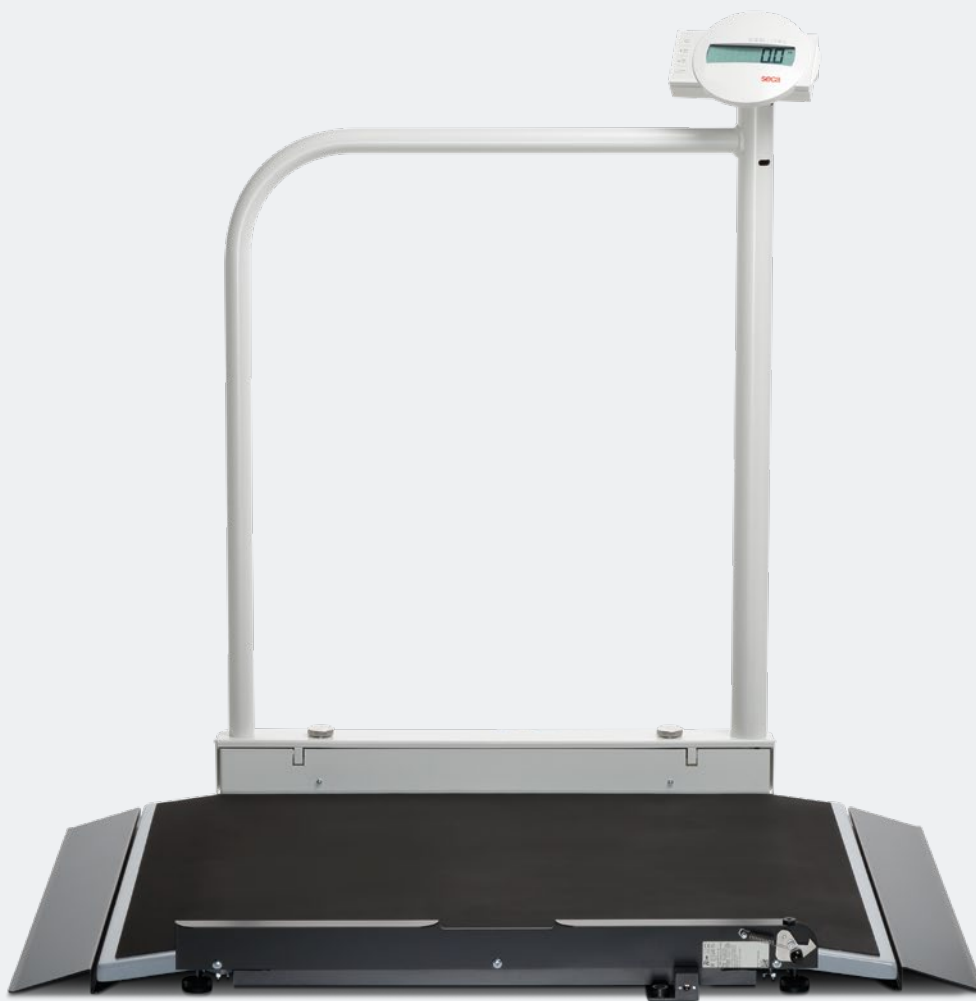

seca 677 r

Wheelchair scale with handrail  
and RS-232 interface



seca 677 r – wheelchair accessible weighing  
made easy

The seca 677 r can be used as a wheelchair scale or as a “seat scale” with a chair, where its tare weight is deducted. The handrail always provides a secure feeling. It is so firmly anchored that even a patient affected by obesity can lean against it without it giving in. The measuring result remains unaffected by this. The non-slip rubber pad and the safety rail ensure additional safety, preventing the wheelchair from rolling over to the side.

Approval class	III
Capacity	300 kg
Graduation	100 g < 200 kg > 200 g
Dimensions (WxHxD)	920 × 1,120 × 1,150 mm
Platform (WxHxD)	800 × 55 × 965 mm
Net weight	38 kg
Power supply	Power adapter
Functions	TARE, Pre-TARE, HOLD, Auto-HOLD, Auto-CLEAR, RESET, Automatic weighing range switch-over, Adjustable damping, BMI
Interfaces	RS-232
Medical Device Regulation (MDR) certified	



Mobile and space-saving

If the patient does not come to the seca 677 r, the seca 677 r just comes to the patient. Folded up and with transport castors, an individual person can move the wheelchair scale to another room since the handrail serves as a handle. In order to save space, the seca 677 r can be folded up and removed when required.



We reserve the right to make modifications.